

Bench scales KERN ECE-N · ECB-N



Plain, mobile, uncomplicated

Features

- **1** High mobility: Battery operation, compact, lightweight construction, recessed grips on the underside and low weight make this balance ideal for use in several locations
- Particularly large weighing plate
- Very fast display: steady weight values within 2 s
- Simple and convenient 2-key operation

Technical data

- KERN ECB-N: Large backlit LCD display, digit height 21 mm
- KERN ECE-N: Large LCD display, digit height 21 mm
- Dimensions weighing surface W×D 320×260 mm
- Weighing plate material
 - KERN ECE-N: plastic
 - KERN ECB-N: stainless steel, see larger picture
- Overall dimensions W×D×H 320×300×60 mm
- Optional battery operation, 6×1.5 V AA, standard, operating time up to 100 h, AUTO-OFF function to preserve the battery
- Net weight
 - KERN ECE-N: approx. 1,6 kg
 - KERN ECB-N: approx. 2,6 kg
- Permissible ambient temperature 5 °C/35 °C


Accessories

- Internal rechargeable battery pack, operating time up to 30 h without backlight, charging time approx. 10 h, KERN PCB-A01
- Rechargeable battery pack external, operating time up to 30 h without backlight, charging time approx. 10 h, KERN KS-A01
- Tare pan made from stainless steel. Ideal for weighing loose small parts as well as fruit, vegetables, etc., KERN RFS-A02
- External universal mains adapter, with universal input and optional input socket adapters for EU, GB, USA, KERN YKA-03
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD



OPTION

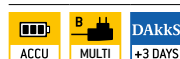

Model	Weighing capacity [Max] kg	Readability [d] g	Reproducibility g	Linearity g	Option	
					DAkkS Calibr. Certificate	
KERN					DAkkS	
					KERN	
ECE 10K-3N	10	5	5	± 20	963-128	
ECE 20K-2N	20	10	10	± 40	963-128	
ECE 50K-2N	50	20	20	± 80	963-128	
<hr/>						
ECB 10K-3N	10	5	5	± 20	963-128	
ECB 20K-2N	20	10	10	± 40	963-128	
ECB 50K-2N	50	20	20	± 80	963-128	

Pictograms

Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)	KERN Communication Protocol (KCP): It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems	Suspended weighing: Load support with hook on the underside of the balance
Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required	GLP/ISO log: The balance displays serial number, user ID, weight, date and time, regardless of a printer connection	Battery operation: Ready for battery operation. The battery type is specified for each device
Easy Touch: Suitable for the connection, data transmission and control through PC or tablet.	GLP/ISO log: With weight, date and time. Only with KERN printers.	Rechargeable battery pack: Rechargeable set
Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	GLP/ISO log: With weight, date and time. Only with KERN printers.	Universal plug-in power supply: with universal input and optional input socket adapters for A) EU, CH, GB; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS
Alibi memory: Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.	Piece counting: Reference quantities selectable. Display can be switched from piece to weight	Plug-in power supply: 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available
Data interface RS-232: To connect the balance to a printer, PC or network	Recipe level A: The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out	Integrated power supply unit: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
RS-485 data interface: To connect the balance to a printer, PC or other peripherals. Suitable for datatransfer over large distances. Network in bus topology is possible	Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display	Weighing principle: Strain gauges: Electrical resistor on an elastic deforming body
USB data interface: To connect the balance to a printer, PC or other peripherals	Totalising level A: The weights of similar items can be added together and the total can be printed out	Weighing principle: Tuning fork: A resonating body is electromagnetically excited, causing it to oscillate
Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals	Percentage determination: Determining the deviation in % from the target value (100 %)	Weighing principle: Electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings
WiFi data interface: To transfer data from the balance to a printer, PC or other peripherals	Weighing units: Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details	Weighing principle: Single cell technology: Advanced version of the force compensation principle with the highest level of precision
Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.	Weighing with tolerance range: (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model	Verification possible: The time required for verification is specified in the pictogram
Analogue interface: to connect a suitable peripheral device for analogue processing of the measurements	Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value	DAKkS calibration possible (DKD): The time required for DAKkS calibration is shown in days in the pictogram
Interface for second balance: For direct connection of a second balance	Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.	Factory calibration (ISO): The time required for Factory calibration is shown in days in the pictogram
Network interface: For connecting the scale to an Ethernet network		Package shipment: The time required for internal shipping preparations is shown in days in the pictogram
		Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg - 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages DE, EN, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Your KERN specialist dealer: